


NEW

RADIANCE® ULTRA TRUCOLOR 32"

Premium Endoscopy Visualization



- Optimized for 1080p endoscopy
- Ultra-wide gamut color calibrated
- Advanced imaging capabilities
- Cleanable splash-proof design
- 10-year scratch-resistant-glass guarantee
- ZeroWire® embedded receiver optional

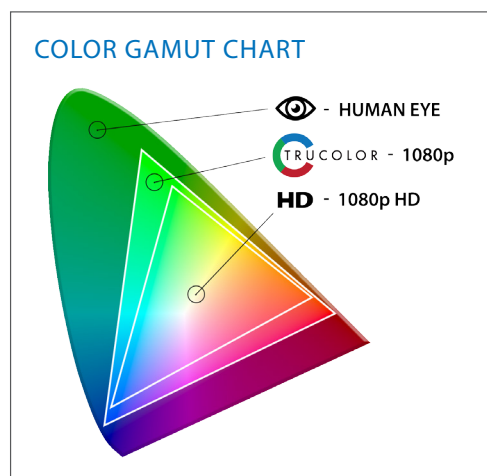

The Radiance® Ultra series leads the industry as a revolutionary display with the aim to transform advanced visualization technology capabilities and features into a clinical solution that supports the drive to improve patient outcomes, improve workflow efficiency, and lower operating costs.

State-of-the-art TruColor technology maximizes a 1080p camera image to provide an endoscopic visualization experience that more accurately resembles anatomy the human eye would visualize during an open surgical procedure. Exclusively engineered over several years by NDS, TruColor technology enhances 1080p video to provide an expanded spectrum of exceptionally vivid colors that is more accurately aligned with the capabilities of human vision. TruColor with Medi-Match™ color calibration assures consistent image quality and accurate color reproduction. The result is a 1080p display with the next-generation color performance of a 4K display.

With a focus to improve workflow efficiency and workplace safety, the Radiance Ultra series is available with an optional built-in ZeroWire receiver. When paired with the ZeroWire *Mobile* battery-powered stand, the combination becomes the world's first and only truly cordless and wireless mobile endoscopic solution (patent pending). To eliminate display scratches caused by IV poles or surgical light heads, the Radiance Ultra series uses scratch-resistant, splash-proof, edge-to-edge glass that includes an industry-exclusive 10-year scratch-resistance guarantee.

RADIANCE® ULTRA TRUCOLOR 32"

Premium Endoscopy Visualization



SPECIFICATIONS

Screen Size	31.5" diagonal
Resolution	1920x1080
Aspect Ratio	16:9
View Angle	178°
Luminance	650 cd/m ² typical @ 6500°K
Backlight Stabilization	NDS Intelli-Guard™
Color Calibration	NDS Medi-Match™
Color Gamut	120% of BT.709
Video Board Configuration	
Primary Board – Digital & Analog	DVI-I and 3G-SDI Input/Output (supports VGA/RGB/SOG)
Secondary Board – Digital	2nd DVI-D and 3G-SDI Inputs/Outputs, Analog Component/VGA
Secondary Board – Analog	2nd DVI-D Input/Output, Analog Component/VGA, S-Video, Composite
Embedded Wireless Receiver	Optional ZeroWire Embedded Receiver
Sealed Glass	Bezel Protected, Edge-to-Edge Glass, Chemically Strengthened
Glass Scratch Warranty	10 Year Scratch-Resistant Guarantee
Display Warranty	3 Years
Dimensions (WxHxD)	30.7" x 20" x 3.4"/780mm x 509mm x 87mm
Weight	29 lbs/12.3 kg (primary board only)
Green Compliance	REACH, RoHS-2, WEEE
Regulatory Compliance	ANSI/AAMI ES60601-1, CAN/CSA C22.2 No. 60601-1, Conflict Minerals, FCC Class B, EN60601-1, EN60601-1-2, CE, MDD 93/42/EEC, Class I Medical Device, Designed to meet IPX6, CCC (Pending)
Regulatory Compliance with ZeroWire Embedded	ANSI/AAMI ES60601-1, CAN/CSA C22.2 No. 60601-1, Conflict Minerals, FCC Class B, CAN ICES-3(B) / NMB-3(B); EN60601-1, EN60601-1-2, MDD 93/42/EEC, EU Class I Medical Device, (R&TTE) Directive 1999/5/EC, this device contains module: SII-SK63101, EN 302 065 V1.2.1, US: Class II Medical Device, FDA cleared, CE0984Φ, FCC ID: UK2-SII-SK63101, IC ID: 6705A-SIISK63101, Designed to meet IPX6

Global Headquarters

www.ndssi.com
5750 Hellyer Ave
San Jose, CA 95138
USA 🇺🇸
+1 408 776 0085
info@ndssi.com

Europe

GSI Group Europe GmbH
Münchner Strasse 2A
82152 Planegg
Germany
+49 89 31 707 100
info.gsig-europe@gsig.com

Asia Pacific

East Square Omori
6-20-14 Minamioi, Shinagawa-ku
Tokyo 140-0013
Japan
+81 3 5753 2466
info@ndssi.jp

NDS Quality System


ISO 13485:2003
FDA Registration #2954921

ORDER INFORMATION

Radiance Ultra TruColor 32"

90R0114	RU TC 32" with Primary Board
90R0112	RU TC 32" with Primary and Secondary Digital Board
90R0113	RU TC 32" with Primary and Secondary Analog Board
90R0115	RU TC 32" with Primary Board and ZeroWire Embedded Receiver

CONFIGURATION

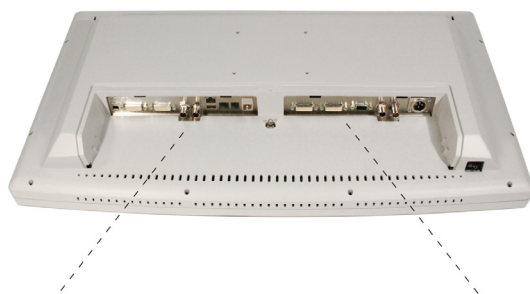

Primary Board

90R0114: Primary Board Only

Secondary Board* (optional)

90R0112: Primary + Digital Board

90R0113: Primary + Analog Board

*Must be installed at factory; not field upgradable.