



The ScaleOR product is a medical-grade video scaling device for use in healthcare environments. It provides the user with the flexibility to choose the appropriate input module per their needs to accommodate either analog or digital signal conversions and scaling. The available input modules are DVI-D, 3G-SDI, RGB/SOG, and S-video/Composite. The video input is scaled, and converted to three simultaneous outputs formats: 3G-SDI, DVI-D and DVI over single fiber (optional). The scaled output format is selectable using a dip-switch setting or by using the on-screen-display menu (OSD). From the OSD menu, users can access advanced features such as RGB gain, brightness, contrast, hue, saturation, zoom, and scaling settings.

➤ Medical-Grade Video Converter/Scaler

➤ Modular Design Supports Both Analog and Digital Video in HD or SD Formats

➤ Multiple Simultaneous HD Video Output Formats


➤ OSD Menu Functionality for Output Format and Other Advanced Features

➤ Universal video format conversion




Medical-Grade Video Scaler

Typical Application Diagram

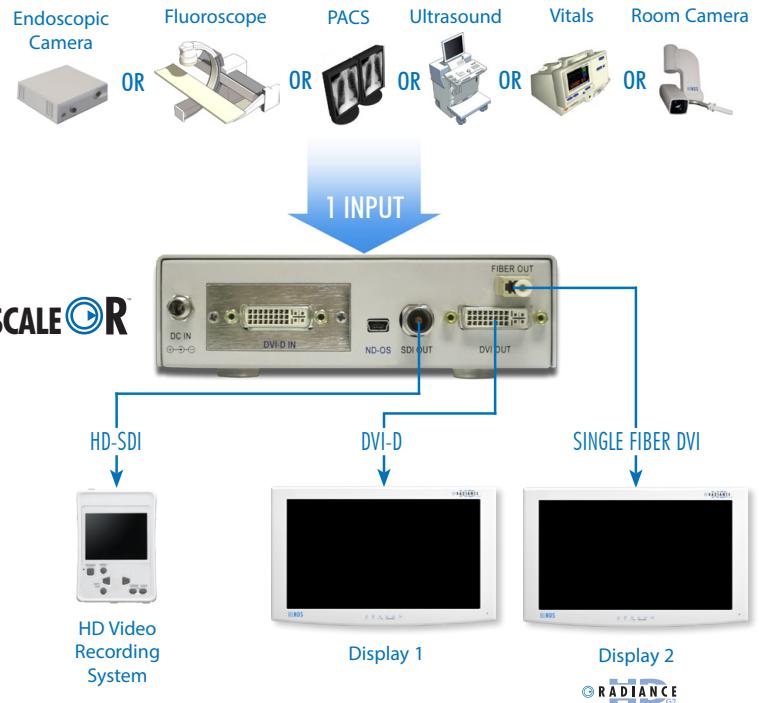

3G-SDI Input Module

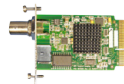

S-Video/Composite Input Module



DVI-D Input Module



VGA/ HD-RGBS/
HD-YP_bP_r/ SOG Input Module

SPECIFICATIONS / FEATURES

Dimensions (W x D x H)	5.80 x 7.5 x 1.67 inches (147.3 x 190.5 x 42.4 mm)
Control Interface	Front Keypad & OSD or serial Interface
Compatible Displays	NDS Radiance®, EndoVue™, and Other HD Displays
Input Signal Detection Indicator	Power Switch Illuminates Steady When Detected
Supports 25 video input formats & 30 graphics resolutions up to 1080p or 1920 x 1200.	
Supports a total of 47 different output resolutions (some on DVI only).	

REGULATORY COMPLIANCE

ANSI/AAMI ES60601-1, CAN/CSA C22.2 No. 60601-1, Conflict Minerals, FCC Part 15, EN60601-1, EN60601-1-2, CE, MDR 2017/745, 2007/47/EC, REACH, RoHS-2, WEEE, Class I Medical Device



Global Headquarters

www.ndssi.com
5750 Hellyer Ave
San Jose, CA 95138
USA 🇺🇸
+1 408 776 0085
info@ndssi.com

Europe

Novanta Europe GmbH
Parkring 57-59
85748 Garching
Germany 🇩🇪
+49 89 31 707 100
info.novanta-europe@novanta.com

Asia Pacific

Novanta Japan
East Square Omori
6-20-14 Minamioi, Shinagawa-ku
Tokyo 140-0013
Japan
+81 3 5753 2466
info@ndssi.com

Oriental Media Center
Suite 2302, Tower C
No.4, Guang Hua Road
Chao Yang District
Beijing, 100026
China
+86 10 8559 7859
info@ndssi.com

ISO 13485
FDA Registration # 2954921